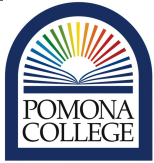


concordance=TRUE

	Environmental Analysis Teaching and Research Laboratory	Date: 8/10/2016	Number: 09 v.01
	Standard Operating Procedure	Title: Using Everyday Lab Equipment	
	Approved By: TBD	Revision Date: 8/10/16	

## Using Everyday Lab Equipment

### 1. Scope and Application

**1.1** The scope of this SOP is to train researchers to learn and follow the basic protocol of using, cleaning, and operating everyday lab equipment, which includes, the micropipette, glassware equipment, measuring balances, the Vortex, and other marginal lab tools that will be further specified.

**1.2** The applications of this SOP are for any basic lab setting for students conducting experiments and/or learning new lab procedures.

### 2. Summary of Method

#### 2.1

#### Contents

<b>1</b>	<b>Scope and Application</b>	<b>1</b>
<b>2</b>	<b>Summary of Method</b>	<b>1</b>
<b>3</b>	<b>Acknowledgements</b>	<b>3</b>
<b>4</b>	<b>Definitions</b>	<b>3</b>
<b>5</b>	<b>Health and Safety</b>	<b>3</b>
<b>6</b>	<b>Personnel &amp; Training Responsibilities</b>	<b>3</b>
<b>7</b>	<b>Apparati</b>	<b>3</b>
<b>8</b>	<b>What is a Micropipette?</b>	<b>3</b>
<b>9</b>	<b>How to Use a Micropipette</b>	<b>3</b>
<b>10</b>	<b>Things to Keep in Mind</b>	<b>4</b>
<b>11</b>	<b>What is a Vortex?</b>	<b>5</b>

12 How and When to Use a Vortex	5
13 QC/QA Criteria	5
14 References	5

### 3. Acknowledgements

### 4. Definitions

#### 4.1 Term1:

### 5. Health and Safety

#### 5.1 Describe the risk...

### Safety and Personnel Protective Equipment

### 6. Personnel & Training Responsibilities

**6.1** Researchers training is required before this the procedures in this method can be used...

**6.2** Researchers using this SOP should be trained for the following SOPs:

- SOP01 Laboratory Safety

### 7. Apparati

#### 7.1 Pipetman (ranging from 2-1000microliters)

### 8. What is a Micropipette?

Micropipettes are used to transfer small amounts (j 1 ml) of liquids. The scales on micropipettes are in microliters (1000 l = 1 ml). These are very expensive, delicate instruments so be very careful during use. The sizes of the micropipettes can be identified via the number on the round button on the plunger; this value is the maximum volume in microliters that can be transferred with that size pipette. The following is an illustration of a micropipette:

### 9. How to Use a Micropipette

1. NEVER exceed the upper or lower limits of these pipettes. The limits are: P10: 1.0 - 10.0 l; P20: 2.0 - 20.0 l; P200: 20 - 200 l; P1000: 200 - 1000 l Look at the front face of the pipet and you will see a window with three digits inside. The diagram below shows the MAXIMUM value that can or should be dialed in on each size pipet. To exceed these values will put the pipet out of calibration. Beside each "window" below is the numbers place it represents. (figure 2)

- **RULE OF THUMB:** Always select the SMALLEST size pipette that will handle the volume you wish to move to achieve the greatest accuracy. Accuracy decreases as you use unnecessarily large pipets for small volumes.
2. Set the desired volume by turning the centrally located rings clockwise to increase volume or counterclockwise to decrease volume.
  3. Load a sterile tip. Use the appropriate tips for each micropipette as designated by the labels of the box, and by confirming with the lab instructor.
  4. 5. Load the sample.
    - The plunger will stop at two different positions when it is depressed. Push the plunger down slowly to the point of first resistance: this is the load volume. Because this first stopping point is dependent on the volume that is being transferred, the distance you have to push the plunger to reach the point of initial resistance will change depending on the volume being pipetted.
    - While holding the plunger at the load volume set point, put the tip into the solution so that it is immersed just enough to cover the end (3-4 mm), not as deep as possible.
    - Slowly release the plunger to draw up the liquid making sure to keep the tip immersed. **NOTE:** If the solution you are pipetting is viscous, allow the pipet tip to fill to final volume before removing it from solution to avoid the presence of bubbles in the plastic tip, which will result in an inaccurate volume.
    - Visually inspect the load to make sure it is correct - there should be no air space in the distal end tip.
  5. Deliver the sample. The second stopping point can be found when the plunger is depressed beyond the initial resistance until it is in contact with the body of the pipette. This second stopping point is used for the complete discharging of solutions from the plastic tip. You should not reach this second stop when drawing liquid into the pipette, only when expelling the last drop.
  6. Discharge the tip. While holding the tip over an appropriate waste receptacle, press the discharge slider on the back of the grip.

## 10. Things to Keep in Mind

- Do not hesitate to ask the lab instructor questions about the pipetting protocol and measures to avoid contamination
- Never point a pipette up. This may cause liquid to run down into the pipette destroying it
- When withdrawing liquids with the pipette, always release the plunger slowly. This prevents liquid from rushing into the end of the pipette and clogging it up. This is especially important with large volume pipettes

- Be sure you use the proper size tip for each pipette.
- Always use a new tip for each different liquid.
- Use the correct pipette for the volume that is to be dispensed. Going below or above the range of the micropipette may damage the instrument.

## 10.1 Vortex Mixer

## 11. What is a Vortex?

A vortex is a relatively simple lab device used primarily to mix small vials of liquid. It consists of an electric motor with a vertical drive shaft attached to a cupped rubber piece mounted slightly off-center. As the motor runs the rubber piece oscillates rapidly in a circular motion. When a test tube or other appropriate container is pressed into the rubber cup (or touched to its edge) the motion is transmitted to the liquid inside and a vortex is created.

## 12. How and When to Use a Vortex

1. Turn on Vortex and either put it in "On" or "Auto" Mode
2. If the Vortex is set to "On" mode, turn the dial to an appropriate speed from 1 to 10. The Vortex will begin spinning at that speed. Then, place the test tube or whatever appropriate container that is holding the sample, gently pressing it against the rubber piece to initiate movement.
3. If it is set on "Auto" mode, the rubber piece will only move when it makes contact with the sample container. This mode is often used because there is more control in terms of spinning the material.

## 12.1 Glassware

## 13. QC/QA Criteria

## 14. References

- 14.1 APHA, AWWA, WEF. (2012) Standard Methods for examination of water and wastewater. 22nd American Public Health Association (Eds.). Washington. 1360 pp. (2014).

Figure 1: The Pipetman

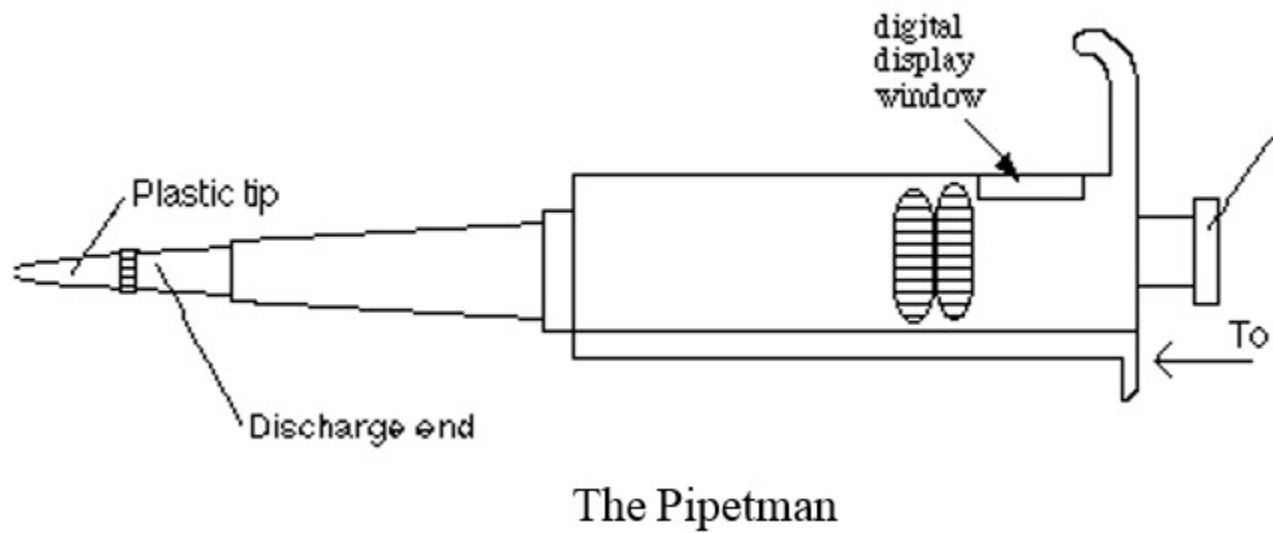


Figure 2: Pipetman Window

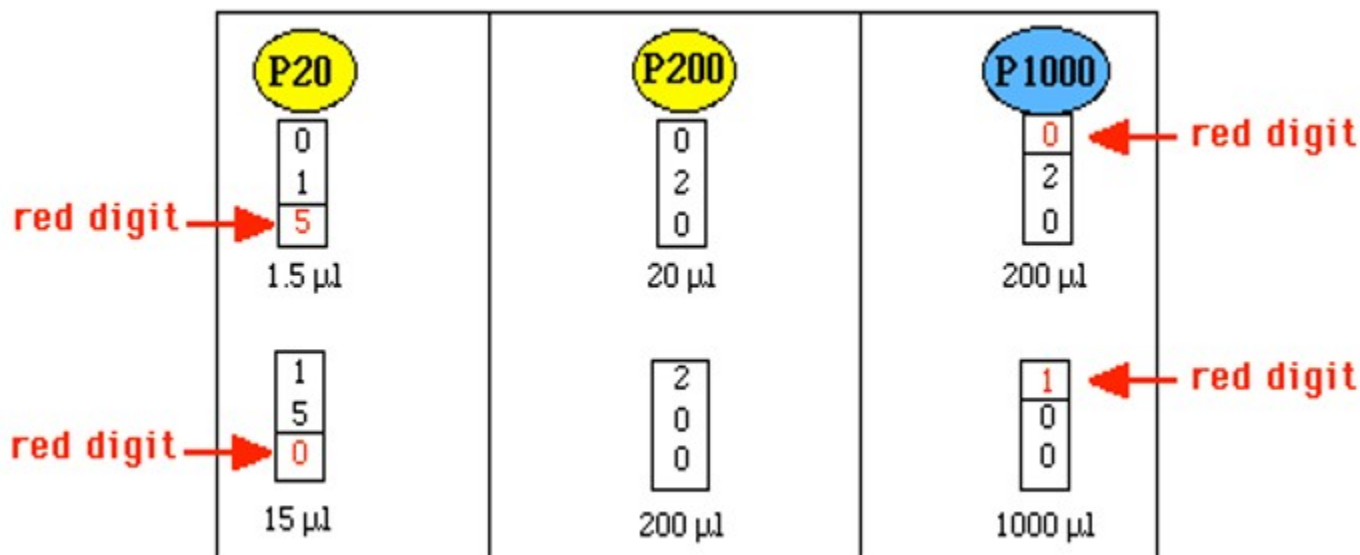


Figure 3: Loading Sample

