

## Dimethylphenol Method

Method 10206<sup>1</sup>0.23 to 13.50 mg/L NO<sub>3</sub><sup>-</sup>-N or 1.00 to 60.00 mg/L NO<sub>3</sub><sup>-</sup> (LR)

TNTplus™ 835

**Scope and application:** For wastewater, drinking water, surface water and process water.<sup>1</sup> USEPA approved for water and wastewater analysis, 40 CFR part 136.**Test preparation****Instrument-specific information**

[Table 1](#) shows all of the instruments that have the program for this test. The table also shows the adapter and light shield requirements for the applicable instruments that can use TNTplus vials.

To use the table, select an instrument, then read across to find the applicable information for this test.

**Table 1 Instrument-specific information for TNTplus vials**

Instrument	Adapters	Light shield
DR 6000, DR 5000	—	—
DR 3900	—	LZV849
DR 3800, DR 2800	—	LZV646
DR 1900	9609900 or 9609800 (A)	—

**Before starting**

DR 3900, DR 3800, DR 2800: Install the light shield in Cell Compartment #2 before this test is started.

Review the safety information and the expiration date on the package.

The recommended sample pH is 3–10.

The sample temperature must be 20–23 °C (68–73 °F) for accurate results.

The recommended temperature for reagent storage is 15–25 °C (59–77 °F).

DR 1900: Go to All Programs>LCK or TNTplus Methods>Options to select the TNTplus number for the test. Other instruments automatically select the method from the barcode on the vial.

Review the Safety Data Sheets (MSDS/SDS) for the chemicals that are used. Use the recommended personal protective equipment.

Dispose of reacted solutions according to local, state and federal regulations. Refer to the Safety Data Sheets for disposal information for unused reagents. Refer to the environmental, health and safety staff for your facility and/or local regulatory agencies for further disposal information.

**Items to collect**


Description	Quantity
Nitrate LR TNTplus Reagent Set	1
Pipet, adjustable volume, 0.2–1.0 mL	1
Pipet tips, for 0.2–1.0 mL pipet	1

Refer to [Consumables and replacement items](#) on page 4 for order information.


## Sample collection and storage

- Collect samples in clean glass or plastic bottles.
- Analyze the samples as soon as possible for best results.
- If immediate analysis is not possible, immediately filter and keep the samples at or below 6 °C (43 °F) for a maximum of 48 hours.
- To preserve samples for a maximum of 14 days, adjust the sample pH to 2 or less with concentrated sulfuric acid (approximately 2 mL per liter) and keep at or below 6 °C (43 °F). The test results then include nitrate and nitrite.
- Let the sample temperature increase to room temperature before analysis.
- Before analysis, adjust the pH to 7 with 5 N sodium hydroxide solution.
- Correct the test result for the dilution caused by the volume additions.


## Test procedure

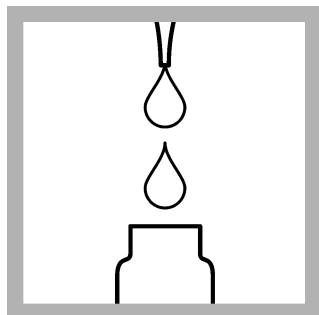

1. Use a pipet to add 1.0 mL of sample to the test vial.

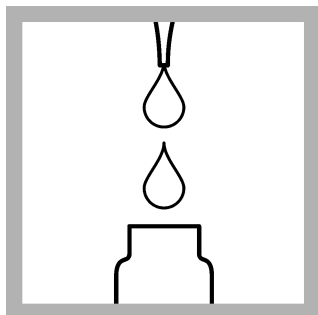


2. Use a pipet to add 0.2 mL of Solution A to the test vial.

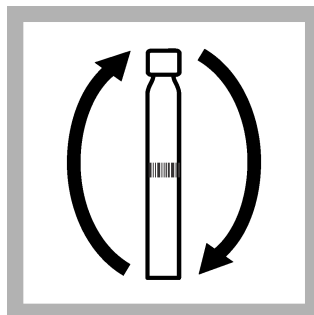


3. Tighten the cap on the vial and invert until completely mixed.




4. Start the reaction time of 15 minutes.

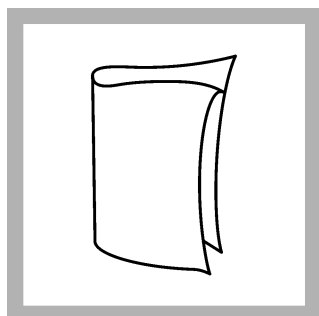

5. When the timer expires, clean the vial.

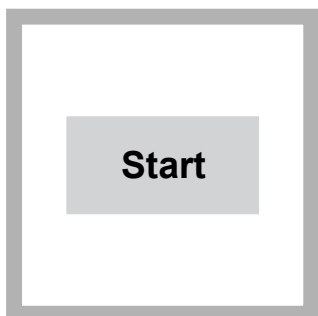

6. DR 1900 only: Select program 835. Refer to [Before starting](#) on page 1.

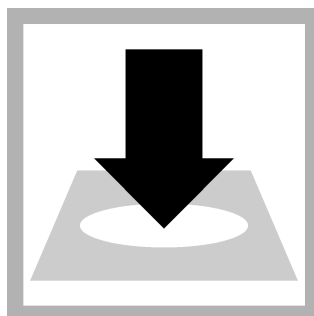

7. Insert the vial into the cell holder. DR 1900 only: Push **READ**. Results show in mg/L  $\text{NO}_3^-$ -N.

## Reagent blank correction

For the best results, measure the reagent blank value for each new lot of reagent. Replace the sample with deionized water in the test procedure to determine the reagent blank value. Subtract the reagent blank value from the sample results automatically with the reagent blank adjust option. Measure the reagent blank value when a new lot of reagent is used.

1. Use deionized water as the sample in the test procedure to measure the reagent blank value.
2. Set the reagent blank function to on. The measured reagent blank value is shown.

3. Accept the blank value. The reagent blank value is then subtracted from all results until the reagent blank function is set to off or a different method is selected.

**Note:** As an alternative, record or enter the reagent blank value at a different time. Push the highlighted reagent blank box and use the keypad to enter the value.

## Sample blanks

If the sample has color or turbidity, measure a sample blank to correct the test result for the interference.

### Items to collect:

- TNTplus 919 sample blank vial
1. Do the test procedure.
  2. Put the sample in the sample blank vial. Fill to the neck of the sample blank vial.
  3. Wipe the sample blank vial clean, then put it into the cell holder. If applicable, the instrument reads the barcode of the sample blank vial and subtracts the value from the initial test result.

## Interferences

Table 2 shows that the ions were individually examined to the given concentrations and do not cause interference. No cumulative effects or influences of other ions were found.

The cumulative effects and influence of other ions have not been found. High loads of oxidizable organic substances (COD) cause the reagent to change color and to give high-bias results. The test can thus only be used for wastewater analyses if the COD is less than 500 mg/L. Verify measurement results with sample dilutions or standard additions.

Nitrite concentrations of more than 2.0 mg/L interfere (high-bias results). Add 50 mg of sulfamic acid (amidosulfonic acid) to 5.0 mL of sample, dissolve and wait for 10 minutes. Analyze the prepared sample as described in the procedure above.

**Table 2 Interfering substances**

Interfering substance	Interference level
Cl <sup>-</sup>	500 mg/L
K <sup>+</sup>	500 mg/L
Na <sup>+</sup>	500 mg/L
Ca <sup>2+</sup>	50 mg/L
Cd <sup>2+</sup>	50 mg/L
Cu <sup>2+</sup>	50 mg/L
Fe <sup>3+</sup>	50 mg/L
Ni <sup>2+</sup>	50 mg/L
Pb <sup>2+</sup>	50 mg/L
Sn <sup>2+</sup>	50 mg/L
Zn <sup>2+</sup>	50 mg/L
Cr <sup>6+</sup>	5 mg/L
NO <sub>2</sub> <sup>-</sup>	2 mg/L
Ag <sup>+</sup>	100 mg/L
Co <sup>2+</sup>	10 mg/L
Fe <sup>2+</sup>	10 mg/L

## Accuracy check

### Standard solution method

Use the standard solution method to validate the test procedure, the reagents and the instrument.

Items to collect:

- Nitrate-Nitrogen Standard Solution, 10.0-mg/L  $\text{NO}_3^-$ -N or Wastewater Influent Standard Solution, Mixed Parameter

1. Use the test procedure to measure the concentration of the standard solution.
2. Compare the expected result to the actual result.

**Note:** The factory calibration can be adjusted slightly with the standard adjust option so that the instrument shows the expected value of the standard solution. The adjusted calibration is then used for all test results. This adjustment can increase the test accuracy when there are small variations in the reagents or instruments.

### Summary of Method

Nitrate ions in solutions that contain sulfuric and phosphoric acids react with 2,6-dimethylphenol to form 4-nitro-2,6-dimethylphenol. The measurement wavelength is 345 nm.

### Consumables and replacement items

#### Required reagents

Description	Quantity/Test	Unit	Item no.
Nitrate LR TNTplus Reagent Set	1	25/pkg	TNT835

#### Required apparatus

Description	Quantity/test	Unit	Item no.
Pipet, adjustable volume, 0.2–1.0 mL	1	each	BBP078
Pipet tips, for 0.2–1.0 mL pipet	2	100/pkg	BBP079
Light shield, DR 3800, DR 2800, DR 2700	1	each	LZV646
Light shield, DR 3900	1	each	LZV849

#### Recommended standards

Description	Unit	Item no.
Nitrate Nitrogen Standard Solution, 10.0-mg/L $\text{NO}_3^-$ -N	500 mL	30749
Nitrate Nitrogen Standard Solution 1000-mg/L $\text{NO}_3^-$ -N	500 mL	1279249
Wastewater Influent Standard Solution, Mixed Parameter, for $\text{NH}_3$ -N, $\text{NO}_3^-$ -N, $\text{PO}_4$ , COD, $\text{SO}_4$ , TOC	500 mL	2833149

#### Optional reagents and apparatus

Description	Unit	Item no.
Filter membrane, 0.45-micron, 25-mm	100/pkg	2514101
Sampling bottle with cap, low density polyethylene, 500-mL	12/pkg	2087079
Sodium Hydroxide Standard Solution, 5.0 N	100 mL MDB	245032
Sulfamic Acid, 454 g	each	234401
Sulfuric Acid, concentrated, ACS	500 mL	97949

---

**Optional reagents and apparatus (continued)**

Description	Unit	Item no.
Test tube rack, polyethylene, for 13-mm OD vials, 90 holes	each	2497900
Water, deionized	4 L	27256



**FOR TECHNICAL ASSISTANCE, PRICE INFORMATION AND ORDERING:**

In the U.S.A. – Call toll-free 800-227-4224

Outside the U.S.A. – Contact the HACH office or distributor serving you.

On the Worldwide Web – [www.hach.com](http://www.hach.com); E-mail – [techhelp@hach.com](mailto:techhelp@hach.com)

**HACH COMPANY**

WORLD HEADQUARTERS

Telephone: (970) 669-3050

FAX: (970) 669-2932